

# Instructions for use



MILLIPORE

## Gel extraction kit Ultrafree<sup>®</sup> DA

DNA Extraction from Agarose Gels, Range: 100 –10.000 bp DNA

AE86.1

**Also suitable for isolation of  
DNA from polyacrylamide gels.  
Use same protocol.**

### Introduction

This centrifugal filter device is designed to extract 100 to 10,000 bp DNA from agarose gel slices in one 10 minute spin. DNA prepared with this device requires no further purification for use in cloning and radioisotopic or fluorescent DNA sequencing. Due to the high resolving power of agarose gel electrophoresis, the small and large non-specific amplification products that frequently interfere with cloning and sequencing after PCR\* are completely removed from your product.

\* PCR (polymerase chain reaction) is covered by U.S. and other equivalent patents issued to Hoffman-La Roche, Inc.

### Contents of Package

### Quantity/pk (preassembled)


Gel Nebulizer, Ultrafree-MC (0.45 µm Durapore)	50
and Microcentrifuge filtrate vials	50

### Additional Materials Required

- TAE electrophoresis buffer (40 mM Tris-acetate, pH 8.0, 0.1 mM Na<sup>2</sup>EDTA).  
**Note:** In principle, the Ultrafree-DA gel extraction kit works with all standard buffer systems like TBE or TAE buffers. However, since gel extractions are generally performed prior to enzymatic processes done with the eluted DNA, we recommend modified TAE (Roth TAE *light*, Art. No. 0122.1) rather than TBE for the following reasons:  
(1) TBE buffer strongly inhibits DNA sequencing reactions while modified TAE buffer does not. (2) Modified TAE has reduced contents of Na<sub>2</sub>-EDTA while regular TAE has 1.0 mM Na<sub>2</sub>-EDTA. This EDTA level will not interfere with the magnesium concentration in sequencing reactions and other downstream enzymatic treatments, many that are dependent on magnesium.
- Agarose for genetic engineering (e.g. Carl Roth Agarose GTQ, Art. No. 6352.2)  
**Note:** Low melting point agarose is not compatible with this protocol.

## Procedure to Use Ultrafree-DA

1. Electrophorese 30 µl of PCR product or other DNA through a < 1.25% ordinary agarose gel, prepared in TAE buffer with ethidium bromide (0.5 µg/mL).
2. Locate band of interest with the long-wavelength ultraviolet lamp or transilluminator. With a razor blade, cut out the slice of agarose (<100 µl at 100 mg) containing the band of interest. The Ultrafree-DA is pre-assembled as follows:
3. Place gel slice into Gel Nebulizer and seal device with the cap attached to vial.
4. Spin at 5,000 x g for 10 minutes. Centrifugation forces the agarose through the Gel Nebulizer, converting it to a fine slurry that is captured by Ultrafree-MC. Extruded DNA in electrophoresis buffer passes through the microporous membrane in Ultrafree-MC and collects in the filtrate vial.
5. Discard the Ultrafree-MC and Gel Nebulizer. You can sequence or clone the DNA in the filtrate without further purification.

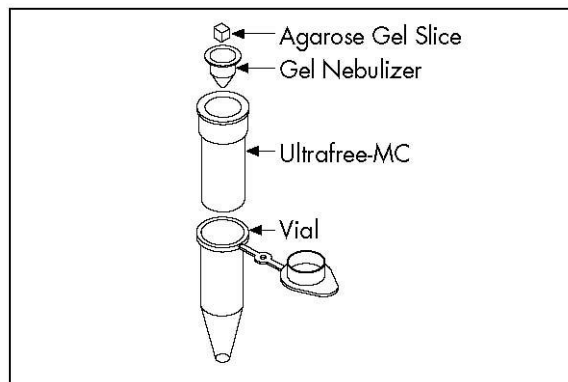

## Technical Assistance

For more information, contact Carl Roth GmbH + Co. KG (Karlsruhe, Germany, 0049 – 721 5606-0) or the Millipore office nearest you.

In the U.S., call **1-800-MILLIPORE** (1-800-645-5476). Outside the U.S., see your Millipore laboratory catalogue for the phone number of the office nearest you. You can reach us by e-mail at [tech\\_service@millipore.com](mailto:tech_service@millipore.com) or visit our web site ([www.millipore.com](http://www.millipore.com)). Millipore Corporation is pleased to provide internet access to Material Safety Data Sheets (MSDS) for its products that contain hazardous materials.

To obtain any MSDS documents that may be associated with this product, go to the MSDS page of the Carl Roth website ([www.carlroth.com](http://www.carlroth.com)) or the Millipore website ([www.millipore.com/msds.nsf/home](http://www.millipore.com/msds.nsf/home)).

## Standard Warranty

**Millipore Corporation** ("Millipore") warrants its products will meet their applicable published specifications when used in accordance with their applicable instructions for a period of one year from shipment of the products. **MILLIPORE MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED. THERE IS NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.** The warranty provided herein and the data, specifications and descriptions of Millipore products appearing in Millipore's published catalogues and product literature may not be altered except by express written agreement signed by an officer of Millipore. Representations, oral or written, which are inconsistent with this warranty or such publications are not authorized and if given, should not be relied upon. In the event of a breach of the foregoing warranty, Millipore's sole obligation shall be to repair or replace, at its option, the applicable product or part thereof, provided the customer notifies Millipore promptly of any such breach. If after exercising reasonable efforts, Millipore is unable to repair or replace the product or part, then Millipore shall refund to the customer all monies paid for such applicable product or part. **MILLIPORE SHALL NOT BE LIABLE FOR CONSEQUENTIAL, INCIDENTAL, SPECIAL OR ANY OTHER INDIRECT DAMAGES RESULTING FROM ECONOMIC LOSS OR PROPERTY DAMAGE SUSTAINED BY ANY CUSTOMER FROM THE USE OF ITS PRODUCTS.** Millipore is a registered trademark of Millipore Corp. or an affiliated company. Ultrafree is a registered trademark of Millipore Corp. or an affiliated company. SeaKem is a registered trademark of FMC BioProducts. 99467 Rev.C, 01/01

**Gel extraction kit Ultrafree® DA**

**50 units**

**AE86.1**

### Carl Roth GmbH + Co. KG

Schoemperlenstraße 3-5 • 76185 Karlsruhe  
P.O. Box 100121 • 76231 Karlsruhe  
Phone: +49 (0) 721/ 5606-0  
Fax: +49 (0) 721/ 5606-149  
info@carlroth.com • www.carlroth.com

sse 06/2021

Die Firma ist eine Kommanditgesellschaft mit Sitz in Karlsruhe, Reg. Gericht Mannheim HRA 100055. Persönlich haftende Gesellschafterin ist die Roth Chemie GmbH mit Sitz in Karlsruhe, Reg. Gericht Mannheim HRB 100428.  
Geschäftsführer: André Houdelet