| | | |
|------------------|------------------|-------------------|
| C-MAG HS4 | C-MAG HS7 | C-MAG HS10 |
| C-MAG HP4 | C-MAG HP7 | C-MAG HP10 |
| C-MAG MS4 | C-MAG MS4 | C-MAG MS10 |

| | | |
|----------------------------------|-----------|----|
| <i>OPERATING INSTRUCTIONS</i> | <i>EN</i> | 3 |
| <i>BETRIEBSANLEITUNG</i> | <i>DE</i> | 10 |
| <i>MODE D'EMPLOI</i> | <i>FR</i> | 17 |
| <i>INDICACIONES DE SEGURIDAD</i> | <i>ES</i> | 24 |



CMAG 032006

Reg.-No. 4343-01

IKA®
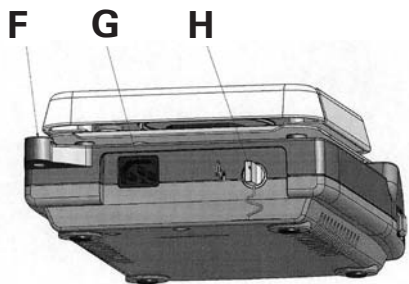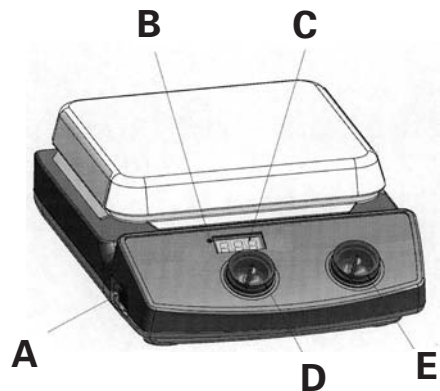
| Item Designation | Pos. Bezeichnung | Pos. Désignation |
|----------------------------|------------------------------|-------------------------------------|
| A Switch | A Geräteschalter | A Commutateur |
| B LED | B LED | B DEL |
| C Control light | C Display | C Affichage |
| D Rotating knob | D Bedienknopf | D Bouton rotatif |
| E Rotating knob | E Bedienknopf | E Bouton rotatif |
| F Threaded support bore | F Stativgewindebohrung | F Alésage fileté du statif |
| G Mains socket | G Netzbuchse | G Prise secteur |
| H Contact-thermometer-jack | H Kontakt-thermometer-Buchse | H douille de thermomètre de contact |

CE-DECLARATION OF CONFORMITY

EN

We declare under our sole responsibility that this product corresponds to the regulations 73/023/EG and 89/336/EG and conforms with the standards or standardized documents EN 61010-1 and EN 61326-1.

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Fig. 1

Contents

| | Page |
|--------------------------------|-------|
| CE - Declaration of conformity | 2 |
| Safety instructions | 3 |
| Unpack | 4 |
| Correct use | 4 |
| Commissioning | 5 |
| Error codes | 6 |
| Accessories | 6 |
| Maintenance | 6 |
| Warranty | 7 |
| Technical data | 8 |
| List of spare parts | 9 |
| Spare parts diagram | 26/27 |

Safety instructions

To your protection


- **Read the operating instructions in full before starting up and follow the safety instructions.**
- Keep the operating instructions in a place where they can be accessed by everyone.
- Ensure that only trained staff work with the appliance.
- Follow the safety instructions, guidelines, occupational health and safety and accident prevention regulations.
- Socket must be earthed (protective ground contact).
- **Attention - Magnetism!** Effects of the magnetic field have to be taken into account (e.g. data cardiac, carriers pacemakers...).
- **Risk of burns!** The heating plate can reach temperatures in excess of 500 °C. Pay attention the residual heat after switching off.
Please make sure that the mains cable does not contact the heating plate
- Wear your personal protective equipment in accordance with the

hazard category of the medium to be processed. Otherwise there is a risk of:

- splashing liquids
- projectile parts
- release any toxic or combustible gases.
- Set up the appliance in a spacious area on an even, stable, clean, non-slip, dry and fireproof surface.
- The feet of the appliance must be clean and undamaged.
- Position the knob at the left stop before starting up. Gradually increase the speed.
- Reduce the speed if
 - the medium splashes out of the vessel because the speed is too high
 - the appliance is not running smoothly
 - the container moves on the set-up surface.
- Never heat up any medium with fire point higher (min. 50°C) than the adjusted temperature
- When using PTFE-coated magnetic bars, the following has to be noted: *Chemical reactions of PTFE occur in contact with molten or dissolved alkaline and alkaline - earth metals, as well as with fine-particled powders of metals of the 2. and 3. group of the periodical system at temperatures above 300-400°C. Only elementary fluorine, chlorine trifluoride und alkaline metals do attack PTFE, halogen hydrocarbons have a reversibly swelling effect.*
Source: Römpps Chemie-Lexikon and „Ullmann“ Bd.19
- Check the appliance and accessories beforehand for damage each time you use them. Do not use damaged components.
- Only replace damaged parts with spare parts identical to the original in function and quality.
- Do not use the device if the ceramic set-up surface is damaged e.g. scratches, splinters or corrosion. A damaged set-up surface could break if used.
- Beware of the risk of
 - flammable materials
 - glass breakage as a result of mechanical shaking power
 - incorrect container size
 - too much medium

- unsafe condition of container
- Only process media that will not react dangerously to the extra energy produced through processing. This also applies to any extra energy produced in other ways, e.g. through light irradiation.
- **Do not** operate the appliance in explosive atmospheres, with hazardous substances or under water.
- A separation from the line is made with the equipment only by pulling net and/or device plug.
- Safe operation is only guaranteed with the accessories described in the "Accessories" chapter.
- Always disconnect the plug before fitting accessories.
- Accessories must be securely attached to the device and cannot come off by themselves. The centre of gravity of the assembly must lie within the the set-up surface.
- The appliance starts up again automatically following a cut in the power supply.
- The appliance may heat up when in use.
- Abrasion of the dispersion equipment or the rotating accessories can get into the medium you are working on.

- To the protection of the equipment** The voltage stated on the nameplate must correspond to the mains voltage.
- Do not coer the device, even partially e.g. with metallic plates or film. This results in overheating.
 - Protect the appliance and accessories from bumps and impacts.

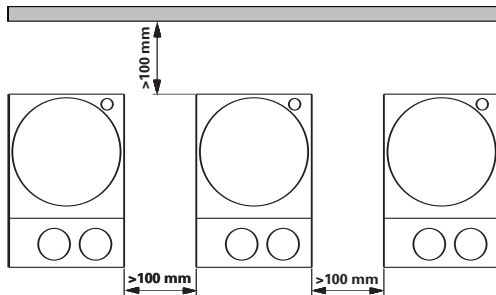- Observe the minimum distances between devices, between the device and the wall and above the assembly (min. 800mm)

Unpack

• Unpack

- Please unpack the device carefully
- In the case of any damage a fact report must be set immediately (post, rail or forwarder)

• Delivery scope

- Heating magnetic stirrer or
Magnetic stirrer or
Heating device
- Mains cable
- Operating instructions

Correct use

• Use

- For mixing and/or heating liquids

• Range of use

- Laboratories
- Chemical industry
- Schools
- Pharmacies

Commissioning

| | MS 4 | MS 7 | MS 10 | HS 4 | HS 7 | HS 10 | HP 4 | HP 7 | HP 10 |
|--|---|------|-------|------|------|-------|------|------|-------|
| Commis- sioning | Put device switch (A) in the OFF position | | | | | | | | |
| | Plug in (G) mains cable | | | | | | | | |
| | Once connected to the power supply the device is in "stand-by" mode | | | | | | | | |
| | The right decimal point on the display (C) is lit | | | | | | | | |
| Stirring | Put device switch (A) in the ON position | | | | | | | | |
| | Any set values are retained when device is switched off and even after the device is disconnected from the power! | | | | | | | | |
| | Set the engine speed with the operating button (E) on the right | | | | | | | | |
| Heating | Put device switch (A) in the ON position | | | | | | | | |
| | Any set values are retained when device is switched off and even after the device is disconnected from the power! | | | | | | | | |
| | Set the target temperature for the heating plate using operating button (D) [for VHP (E)] | | | | | | | | |
| | The set value is indicated on the display (C). If energy is being supplied to the heating plate, the red LED (B) lights up | | | | | | | | |
| | In the stir and stand-by modes, "hot" flashes on the display after the heat has been switched off for as long as the temperature of the set-up surface exceeds 50°C | | | | | | | | |
| Contact- thermometer connection | Put device switch (A) in the OFF position, unplug contact plug | | | | | | | | |
| | Safety contact thermometer acc. to DIN 12878 class 2 connected with jack (H) | | | | | | | | |
| | Put device switch (A) in the ON position | | | | | | | | |
| | Beware the instruction manual of the contact thermometer | | | | | | | | |
| | Important: Display (C) also indicates the target temperature of the heating plate when the contact thermometer is connected | | | | | | | | |

Error codes

Footstep during the enterprise a disturbance up, is indicated these with the devices to **HS** and **HP** by an error message in the display (C). Whenever an error message appears switch the device off to cool down.

| Error code | Cause | Solution |
|------------|----------------------------|---|
| E1 | Inner temperature too high | <ul style="list-style-type: none"> switch off the unit and allow it cool down |
| E6 | Motor seized | <ul style="list-style-type: none"> Stirring bars ungeeignet Reduce the viscosity Foreign bodies in the equipment remove |
| E9 | Safety circuit ETC | <ul style="list-style-type: none"> Plug the contact plug (G) Plug the contact thermometer Connecting cable, plug or contact thermometer defect, exchange |

Have the device repaired if the error is not corrected using the measures described or if another error code is displayed.

Accessories

Stirring organs

stirring bars: Ø 6 mm, lengths to 15mm
 Ø 8 mm, lengths to 50mm
 Ø10mm, length to 80mm

Any other accessories

RSE stirring bar remover
 H15 bath top
 H28 bath top
 H36 holding rod
 H16V support rod

H44
 ETS-D4

cross sleeve
 electrical contact thermometer

Installing the support rod

The support rod is attached using the threaded support bore (F).

- Screw nut M10 on to the support rod as far as the stop
- Screw on the support rod as far as the stop by hand
- Tighten the support rod and nut M10 using a flat wrench (SW17).
- Use bossheads to assemble accessory parts or accessory devices.


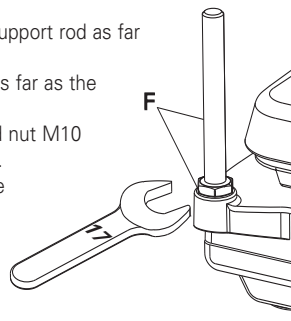
Fig. 3

Maintenance

The appliance is maintenance-free.

Cleaning

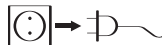Only use cleansing agents which have been recommended by **IKA**.

Use to remove:

| | |
|------------------------|---|
| Dyes | isopropyl alcohol |
| Construction materials | water containing tenside/ isopropyl alcohol |
| Cosmetics | water containing tenside/ isopropyl alcohol |
| Foodstuffs | water containing tenside |
| Fuels | water containing tenside |

- Do not allow moisture to get into the appliance when cleaning
- Wear protective gloves during cleaning the devices.
- Before using another than the recommended method for cleaning or decontamination, the user must ascertain with **IKA** that this method does not destroy the instrument.
- The rubber of the attachments' pads (VG 3.2 and VG 3.3) is a natural product, so talcum can discharge. Clean the rubber pad with glycerin.

Spare parts order

When ordering spare parts, please give:

- Machine type
- Manufacturing number, see type plate
- Item and designation of the spare part, see spare parts list and diagram

Repair

In case of repair the device has to be cleaned and free from any materials which may constitute a health hazard.

If you require servicing, return the appliance in its original packaging. Storage packaging is not sufficient. Please also use suitable transport packaging.

Warranty

In accordance with **IKA** warranty conditions, the warranty period is 24 months. For claims under the warranty please contact your local dealer. You may also send the machine direct to our works, enclosing the delivery invoice and giving reasons for the claim. You will be liable for freight costs.

The warranty does not cover wearing parts, nor does it apply to faults resulting from improper use or insufficient care and maintenance contrary to the instructions in this operating manual.

Technical data

| | | MS 4 | HS 4 | HP 4 | MS 7 | HS 7 | HP 7 | MS 10 | HS 10 | HP 10 | |
|--------------------------------------|------------|-----------------|----------|---------|-----------------|----------|------------|-----------------|----------|---------|--|
| Operating voltage | VAC | | | | | | 230 ±10% | | | | |
| | VAC | | | | | | 120 ±10% | | | | |
| | VAC | | | | | | 100 ±10% | | | | |
| Nominal voltage | VAC | | | | | | 230 / 50Hz | | | | |
| | VAC | | | | | | 120/ 60 Hz | | | | |
| | VAC | | | | | | 100/ 60 Hz | | | | |
| Design frequency | Hz | 50/60 | | | | | | | | | |
| Input power max. at 230 and 120 VAC | W | 30 | 270 | 255 | 30 | 1020 | 1005 | 30 | 1520 | 1505 | |
| | 100 VAC | 30 | 270 | 255 | 30 | 1020 | 1005 | 30 | 1070 | 1055 | |
| Power consumption in stand-by mode | W | 2.5 | | | | | | | | | |
| Perm. duration of operation | % | 100 | | | | | | | | | |
| Perm. ambient temperature | °C | +5 bis +40 | | | | | | | | | |
| Perm. relative humidity | % | 80 | | | | | | | | | |
| Protection type acc. to DIN EN 60529 | | IP 21 | | | | | | | | | |
| Protection class | | I | | | | | | | | | |
| Overvoltage categorie | | II | | | | | | | | | |
| Contamination level | | 2 | | | | | | | | | |
| Operation at a terrestrial altitude | m | max. 2000 | | | | | | | | | |
| Dimensions (W x D x H) | mm | 150 x 260 x 105 | | | 220 x 335 x 105 | | | 300 x 415 x 105 | | | |
| Weight | kg | 3 | | | 5 | | | 6 | | | |
| Motor | | | | | | | | | | | |
| Speed range (infinitely) | 1/min | 100-1500 | 100-1500 | - | 100-1500 | 100-1500 | - | 100-1500 | 100-1500 | - | |
| Speed display | | Scale | Scale | - | Scale | Scale | - | Scale | Scale | - | |
| Power consumption | W | 15 | 15 | - | 15 | 15 | - | 15 | 15 | - | |
| Power output | W | 1,5 | 1,5 | - | 1,5 | 1,5 | - | 1,5 | 1,5 | - | |
| Max. stirring quantity (water) | ltr | 5 | 5 | - | 10 | 10 | - | 15 | 15 | - | |
| Heating plate | | | | | | | | | | | |
| Heating plate dimension | mm | 120x120 | 120x120 | 120x120 | 200x200 | 200x200 | 200x200 | 280x280 | 280x280 | 280x280 | |
| Heating power at 230 and 120 VAC | W | - | 250 | 250 | - | 1000 | 1000 | - | 1500 | 1500 | |
| | at 100 VAC | - | 250 | 250 | - | 1000 | 1000 | - | 1050 | 1050 | |
| Surface temperature min. | °C | - | 50 | 50 | - | 50 | 50 | - | 50 | 50 | |
| Surface temperature max. | °C | - | 500 | 500 | - | 500 | 500 | - | 500 | 500 | |
| Temperature fluctuation | °C | - | ±5 | ±5 | - | ±5 | ±5 | - | ±5 | ±5 | |
| Limit of safety temperature | °C | - | 550 | 550 | - | 550 | 550 | - | 550 | 550 | |

List of spare parts

| Pos. | Designation | MS-C4/C7/C10 | HP-C4/C7/C10 | HS-C4/C7/C10 |
|------|--|--------------|--------------|--------------|
| 7 | Contersunk screw | X | X | X |
| 8 | Connect distributor | X | X | X |
| 12 | Insulating | | X | X |
| 13 | Danger sign | | X | X |
| 17 | Distance holderLP | X | X | X |
| 18 | Clip-Triac fastening | X | X | X |
| 19 | Angel connector | X | X | X |
| 20 | Cover | X | X | X |
| 21 | Mounting plug | X | X | X |
| 31 | Placing plate | X | X | X |
| 32 | Heating foil | | X | X |
| 33 | Insulating | | X | X |
| 35 | Ceramic bushing | | X | X |
| 37 | Holding bracket | X | X | X |
| 38 | Protection sheet | X | X | X |
| 40 | Distance bushing | X | X | X |
| 41 | Distance bushing | X | X | X |
| 45 | Mylar-foil | X | X | X |
| 46 | Color filter foil | | X | X |
| 47 | Front foil | X | X | X |
| 51 | Control knob | X | X | X |
| 53 | Contact plug | | X | X |
| 54 | Magnet | X | | X |
| 59 | Shaded pole asynchron motor | X | | X |
| 65 | Device pedestal | X | X | X |
| 66 | Impulse disk | X | | X |
| 70 | Distance bushing | X | | X |
| 71 | Magnett holder | X | | X |
| 75 | Temperature sensor | | X | X |
| 76 | Pressure sheet | X | X | X |
| 80 | Contact thermometer bushing | | X | X |
| 80 | Cap (Contact-thermometer) | X | | |
| 81 | Connecting cable | X | X | X |
| 100 | Label | X | X | X |
| 2001 | PCB compl. | X | X | X |
| 2002 | PCB accessories consisting of: Poti (motor); Poti (heating); Opto- coupler and LED-display | X | X | X |

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